



Diode Laser System

ICE H8



THE UNPARALLELED STAR FOR
**PROFESSIONAL HAIR
REMOVAL**

Grow Your Business With ICE H8



3 808/760/1064nm
Wavelength Diode Laser

3 Spot Sizes Handles

1 Advanced Cooling
Technology

More and more people desire to have permanent hair reduction.
They want a treatment that is effective, fast and comfortable.
To meet these needs, we developed this laser, so you can:

Benefit from a rapid return on investment

Increase your patient satisfaction

Grow your business

In order to provide you with optimal treatment possibilities, the portable laser ICE H8 comes with:

- Versatile
808nm/808nm+760nm+1064m diode laser
3 spot sizes handpieces
- Advanced cooling technology


The unique features of the Laser ICE H8 enable you to provide your patients with:

- Maximal treatment comfort
- Long lasting results
- Suitable for a variety of skin types

Technology Behind ICE H8

Laser Power 1000w


12*24mm 12*16mm 12*12mm

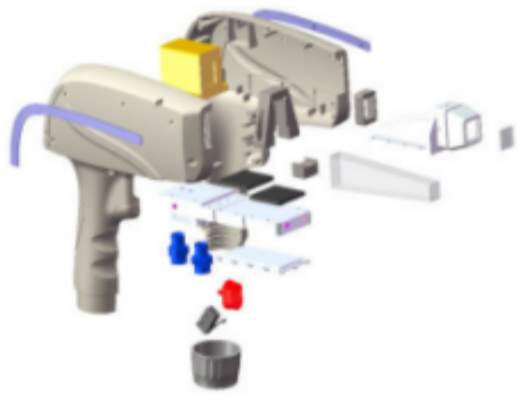

Unique design, high shoot efficiency

Cutting Edge Technology to ensure Painless Hair Removal


Sapphire Touch Cooling
TEC Semiconductor Cooling
Water Cooling

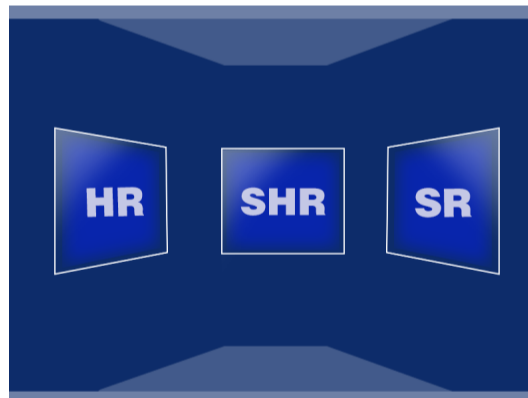

Advanced Graphical User Interface


HR,SHR and SR three functions

Multi-Language optional

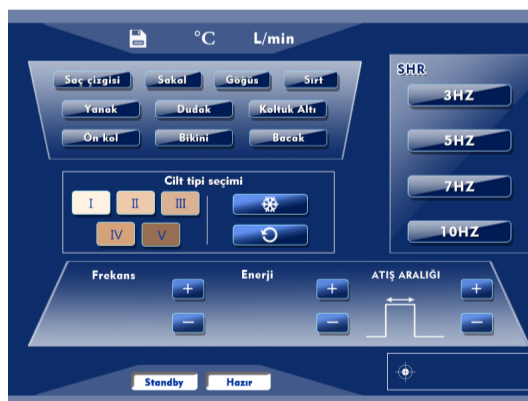


Using the intuitive touch screen, you can select the required mode and programmes.

In each mode (HR or SHR or SR) you can adjust the settings precisely for the skin and hair type and the intensity to obtain the required values for each treatment.

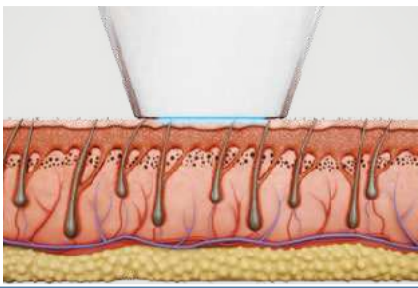

1.Skin is cooled before and during treatment.

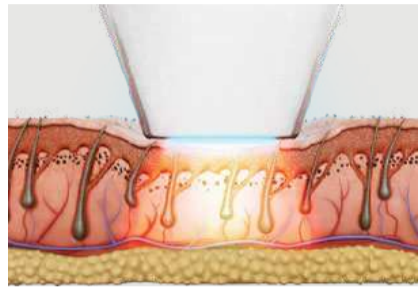

2.Skin is compressed.
A laser beam is applied to the skin and the target is damaged



Cooling System
Water Chiller ,can keep the machine working continuously for 12 hours.



Water Pump
Imported from germany



Water Filter
Effectively filter >15µm imoulties in water



Constant Current Source
Balance current peaks to ensure laser life



Acrylic Water Tank
Avoid corrosion and oxidation of the water tank

Water Line and Circuit Separation
Line optimization design
Protect appliances
Easy maintenance

Meet with You Because of Beauty



Specification

Laser Typ	Diode Laser ICE H8
Wavelength	808nm/808nm+760nm+1064nm
Energy	1-100J/cm2
Treatment head	Sapphire crystal
Pulse Duration	1-300ms
Repetition Rate	1-10 Hz
Interface	10.4
Output power	2000W
Spot Size	12×12mm/12×16mm/12×24mm
Cooling system	Water Chiller cooling +water cooling system+air cooling system+ Sapphire touch cooling
Power Supply	AC 220V 50Hz/ AC 110V 60Hz
Dimension	720*520*600mm
G.W.	43 KG

Spare Parts



Handpiece



Project Glasses



Project Goggles



Foot Switch



Handle Bracket



Filter



Power Cable



Funnel




Key

Warnings and Risks:

Treatment should not be attempted on patients with the following conditions in the treatment area: any active infections or inflammatory skin conditions, Dysplastic nevi, tattoos, active cold sores, open lacerations or abrasions, viral, fungal, or bacterial diseases or a history of post inflammatory hyperpigmentation. The most common immediate responses of laser treatment for hair removal, vascular lesions, pigmented lesions are erythema, edema, perifollicular edema and perivascular edema, vascular bleaching, hyperpigmentation and hypopigmentation or textural changes of the treated areas.

BENEFITS OF LASER HAIR REMOVAL

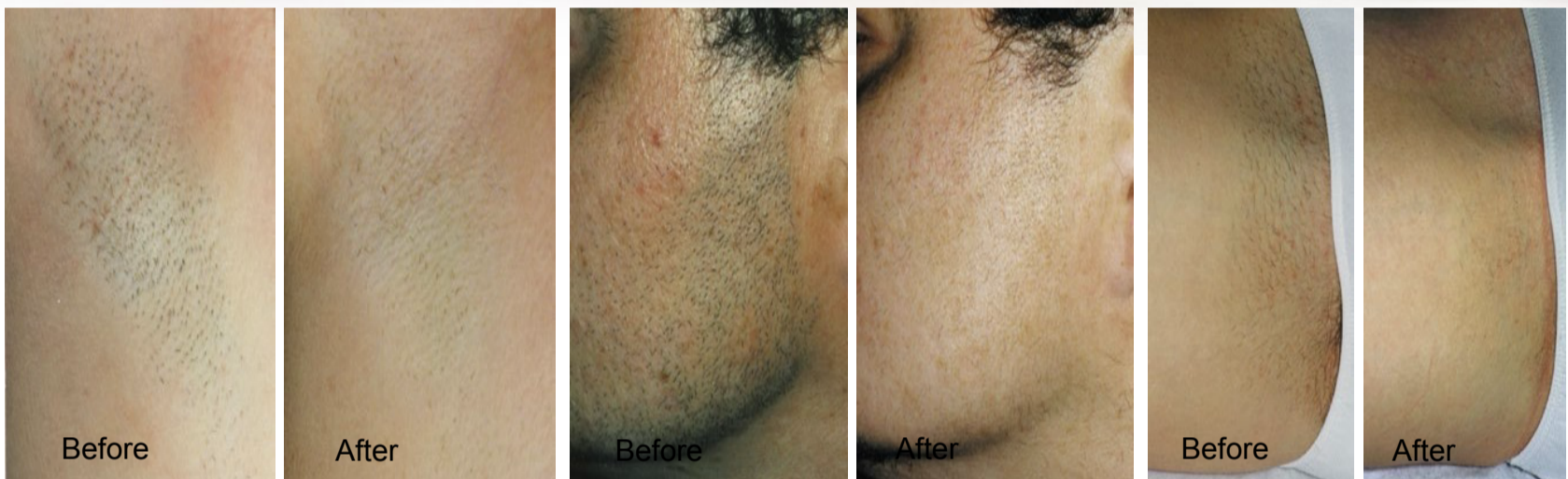

Eight months follow-up post 3 treatments Five months follow-up post 3 treatments Four months follow-up post 3 treatments

“I chose the laser ICE H8 for my practice because it uses an 808 diode and I am very pleased with my choice.

The results are outstanding and the treatment is most comfortable for the patients.”

Adele Quintana, M.D.

“The laser ICE H8 is a fast and easy to use device. It is perfect for treating all body areas. The cooling handpieces enables an effective and comfortable treatment for all of my patients.”

Omar A. Ibrahim, M.D., Ph.D.

TRIANGEL RSD LIMITED

Room 302, Building 3, Weiwei Science and Technology Park, Lekai North Street, Jingxiu District, Baoding City, 071000 China

Email: Yolanda@triangelaser.com

Web: www.triangelaser.com

Whatsapp/Wechat: 0086-15033841027

Europe office for after-sale service:

ADD: Praceta António Gedeão loja 3A

2635- 002 Vale Mourão Portugal

Phone: +351-213540102

Email: geral@iberbraz.com

Trainer and Representative in Europe:

Company: ANABELA PEREIRA AFONSO VAZ

Phone: +351-964482806

Email: anabela_apav@hotmail.com